

# E-Scope II<sup>®</sup>

Improve the way you work

By improving the way you listen

Introducing a new type of stethoscope



- 📻 **Sounds** - The maximum volume of the E-Scope II is approximately 30 times louder than an acoustic stethoscope!
- 📻 **New Specialist Adult Diaphragm** - new stainless steel chest piece included with every E-Scope
- 📻 **Lifetime Technical Support** - call us and a live technician answers the phone
- 📻 **2 minute timer** - for battery conservation
- 📻 **Models for Clinical use, EMS, hearing impaired** - the most complete product line



*E-Scope Hearing Impaired  
Model for ITE hearing aids*

- 📻 **Listen** - to the sounds through headphones
- 📻 **Belt Clip** - swivels for ease of use
- 📻 **Use over hearing aids** - or for comfort







The **E-Scope Electronic Stethoscope** utilizes state of the art electronic technology to amplify heart, breath and Korotkoff sounds so they become much **easier to hear**. The sound is pure and clean. Compared to acoustic stethoscopes, murmurs and faint blows are much more distinguishable (see the below chart).

The E-Scope is **easy to use**. All controls are on the chest piece of the stethoscope. The controls include:

-  Volume Up
-  Volume Down
-  Heart or Breath Sounds switch

Afraid you might leave the E-Scope ON and lose **battery life**? The E-Scope has a built in timer that automatically shuts off the unit after 2 minutes. To reactivate the E-Scope simply press any control button on the head and you're back in business. The E-Scope remembers the last volume setting upon restart. Battery life is estimated at approximately 6 months using an off-the-shelf AAA battery!

**Additional features** include:

-  External output for recording or telemedicine
-  3 types of eartips
-  Softpack carrying case
-  Operator's manual & instructional CD-ROM

**Hearing Impaired?** The E-Scope Belt Model version is made just for you. Simply place the supplied pair of headphones over your in-the-ear hearing aids, or via direct audio input cable for behind-the-ear hearing aids and hear heart/breath sounds with unparalleled clarity.

**EMS - paramedic?** The E-Scope EMS provides exceptional sound quality in a very high noise area.

For more information, uses and updates visit us at [www.cardionics.com](http://www.cardionics.com), or call us at 800-364-5901.



*E-Scope EMS Electronic Stethoscope*



*The 'Specialist' Heads - Pediatric Bell , Adult Bell, Adult Diaphragm*

The 'Specialist' heads consist of a stainless steel bell and diaphragm with a Delrin non-chill ring. They provide a 50% increase in sound sensitivity, a reduction in outside noise and comfortable feel.

### Ordering Information

Cat. No.

**718-7700** - Clinical Model E-Scope with Specialist adult diaphragm, chest piece/earpiece kit and carrying case, 1 year warranty

**718-7740** - Same as cat no. 718-7700 except supplied with 'Specialist' Bell and Specialist Adult diaphragm

**718-7710** - **Hearing Impaired** Belt Model E-Scope with Specialist Adult Diaphragm, chest piece/ear piece kit, carrying case, headphone, 1 year warranty

**718-7750** - Same as cat no. 718-7710 except supplied with 'Specialist' Adult Bell & Specialist Adult Diaphragm

**718-7800** - E-Scope EMS for use in a noisy ambulance  
*Additional Equipment*


**360-7750** - Specialist Adult Bell

**360-7755** - Specialist Pediatric Bell

**718-0405** - External Headphones for second listener or hearing impaired

**718-7030** - Pocket Monitor software for PDA or PC - displays, records & measures duration and frequency of murmurs.

### E-Scope vs Acoustic Stethoscope

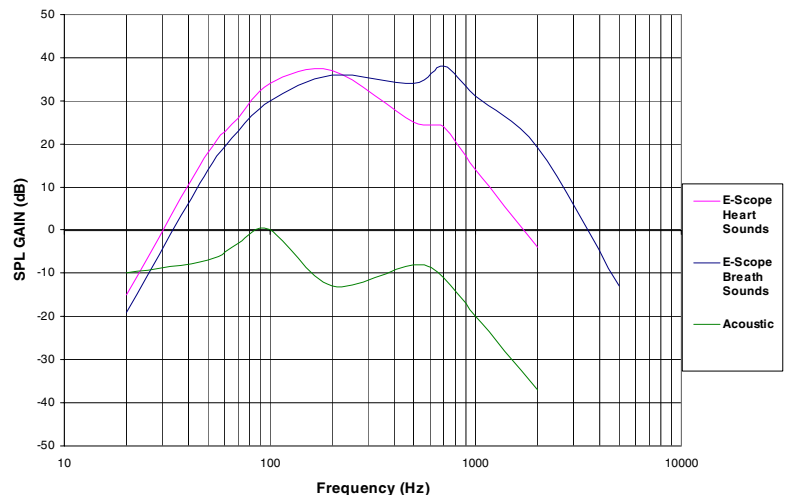

## Cardionics, Inc.

910 Bay Star Blvd.

Telephone : 281-488-5901

Fax: 281-488-3195

Webster, Texas 77598 USA

US and Canada : 1-800-364-5901

Email: keithj@cardionics.com Web site: [www.cardionics.com](http://www.cardionics.com)